





INVACARE **Proflex**

Affordable Pressure Reducing Range

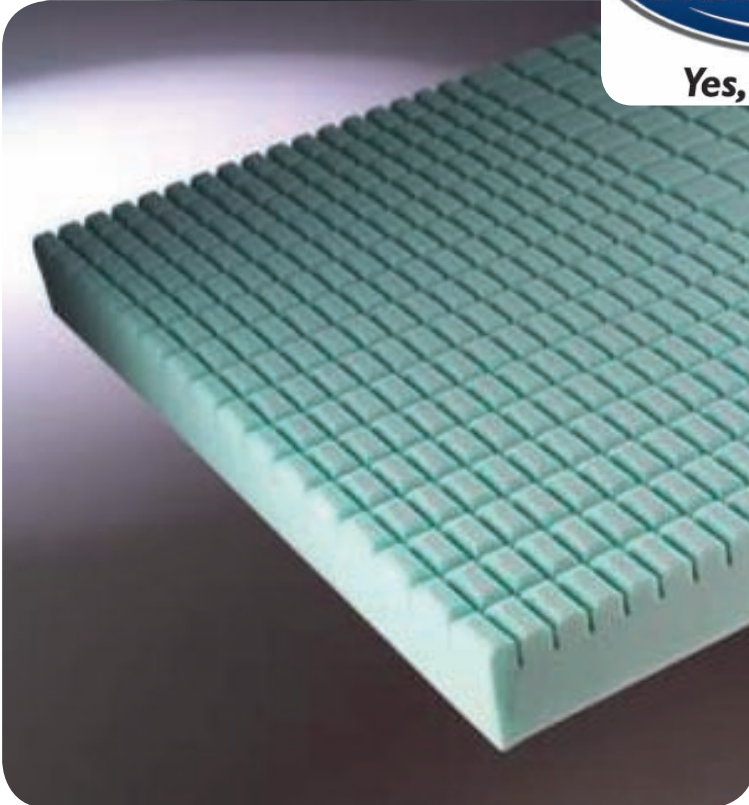

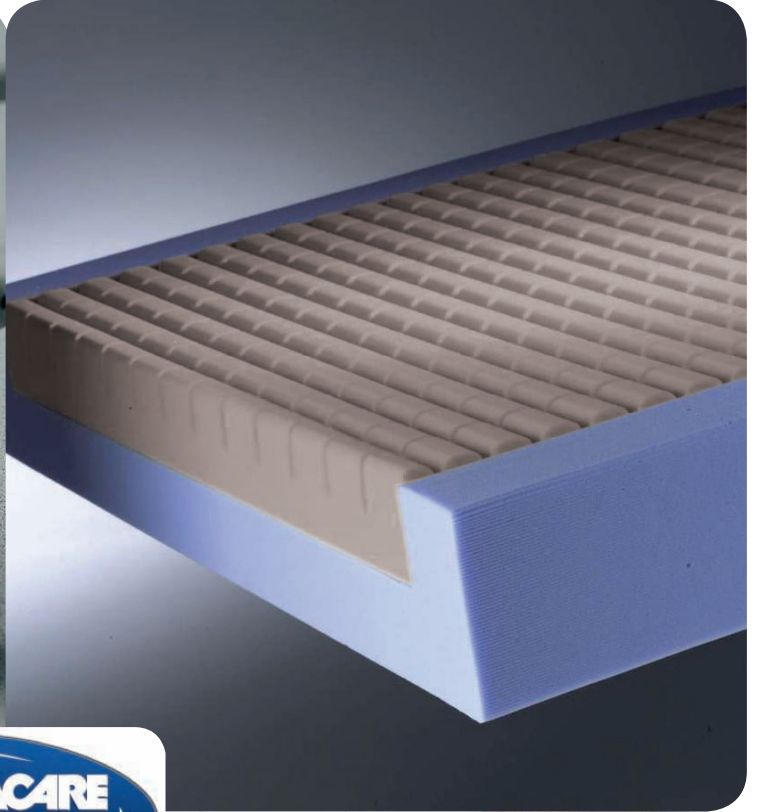

Proflex Wave Visco Cushion

The surface of the Proflex Wave Visco cushion is constructed from moulded visco elastic foam, which is temperature sensitive and conforms excellently to the body, promoting good posture and reducing pressure. The Proflex Visco Wave is suitable for those at ‘High Risk’ of developing pressure ulcers.

Technical Specifications

Risk Category	High Risk*
Weight Limit	20 Stone/127 kg
Warranty	3 Year Foam/ 1 Year Cover
Dimensions	L16 x W16 x D3 ¾” L17 x W17 x D3 ¾” L18 x W16 x D3 ¾” L20 x W18 x D3 ¾” (other sizes available on request)
Replacements	Foam and Cover
Fire Retardancy	BS 7176 (Crib 5)



Features	Benefits
Temperature sensitive and moulded visco elastic foam	Conforms with the body, improving positioning and reducing pressure
Contoured surface	Increases seating surface area to redistribute and reduce pressure
Removable covers	Can be removed for laundering at 80° C or cover replacement.
2 way stretch vapour permeable cover with non slip base	Helps to reduce shear and friction and moisture build up, base reduces slippage

* Essential nursing care is pivotal in pressure ulcer prevention. These products will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



www.invacare.co.uk

www.invacare.ie

Invacare Limited
Pencoed Technology Park
Pencoed
Bridgend
CF35 5HZ
Tel: +44 | 656 776 222
Fax: +44 | 656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Limited
Unit 5
Seatown Business Campus
Seatown Road
SWORDS
County Dublin
Tel: +353 | 810 7084
Fax: +353 | 810 7085
E-mail: ireland@invacare.com
www.invacare.ie

© 2010 Invacare® International Sarl
All rights reserved.

All information quoted is believed to be correct at time of print.
Invacare® reserves the right to alter product specifications without prior consultation.

New Proflex Range - UK
10/2010 - 1539111

