



Invacare® | **Softform® Premier Fully Sealed**

Clinically proven to perform in high risk,
challenging environments



Invacare's fully sealed (FS) pressure reducing mattress meets the needs of challenging environments, such as Mental Health, Secure Units and Dementia Care, where entrapment, ligature risk and deliberate equipment damage such as fire are potential issues. Additionally, the fully sealed nature reduces the risk of items being deliberately concealed inside the mattress cover.

The Fully Sealed (FS) mattress cover is constructed from high specification polyurethane fabric providing multi-stretch, vapour permeable and water proof properties. As a key requirement in specialist areas of healthcare, the cover is manufactured with a Crib 7 rating (BS7177:Crib 7) for fire retardancy.

Effective pressure ulcer prevention

The *Invacare Softform Premier FS* offers the same pressure redistribution benefits as the clinically proven *Softform Premier*. The *Softform Premier* is widely used across the NHS in pressure ulcer prevention programmes. The *Softform Premier FS* combines this clinical efficacy with additional design features that ensure it meets Crib 7 'composite' requirements. **

The multi stretch fabric compliments the internal foam components, conforming effectively to body contours. Combined with the benefit of breathable fabric to enhance microclimate, optimum conditions for pressure ulcer prevention may be achieved.

Specifically designed for high risk areas

The *Invacare FS mattress cover* is made from a highly resistant polyurethane polymer making it extremely resistant to damage caused by harsh chemicals, as well as patient transfer aids. The specialist coating on the FS mattress material provides an impenetrable barrier to viruses and bacteria positively contributing to the challenges of infection control. For areas where the risk of fire is elevated or for areas where it is more difficult for Emergency Services to gain instant access, the Crib 7 FS cover is an ideal solution.

In the absence of zip and fasteners, the FS mattress has fully welded seams. This greatly reduces the risk of intentional penetration or destruction of the internal mattress components.

Product Features




Fully welded seams

The impenetrable fully high frequency welded seams span the perimeter of the mattress ensuring internal components maintain integrity and cannot be accessed.



Each corner has been discretely high frequency welded to produce a neat finish and effective closure to protect internal components from bacteria or fluid ingress.

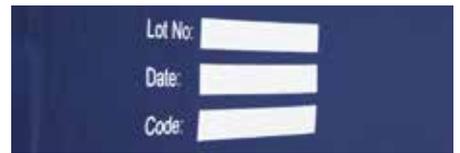

Information Label

Labelling information is screen printed onto the side of the FS cover, eliminating the need for additional material attachments.

Technical data

For more information about this product, including the product's user manual, please visit www.invacare.co.uk

Width	Length	Height	Total product weight	Max. user weight	Max. cleaning temperature	Colour
 85.5-92 cm	 191-200 cm	 15 cm	 13 kg	 247.6 kg (39 stone)	 80°	 Blue

Softform Premier Fully Sealed

A Crib 5 version of the Softform Premier is also available. More information is available on request.

Manufactured to European Standards:
BS EN 597-1 Cig,
BS EN 597-2 Match.

Warranty
1 year

* Essential nursing care is pivotal in pressure ulcer prevention. When selected with a pressure reducing interior, this mattress cover will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value. The assessment tool used in this case is the Waterlow scale.

** Crib 7 'composite' means that the complete mattress components in combination (foam and cover) have been tested and approved by an authorised body (BS7177 CRIB 7, also includes Unintended Use).

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS -
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

Trusted Solutions,
Passionate People

