



Setting New Standards in Bariatric Care

With a safe working load up to 60 Stone (385 kgs), the *Invacare Octave* is a robust yet stylish profiling bed solution that assists in the care and comfort of bariatric clients. With a fully electrical, 4-sectioned profiling mattress support as standard, the **Octave** provides optimal comfort and support to the client, as well as an ideal working environment for the carer. The frameless construction of the mattress platform allows improved side access to the client, and the bed ends can be detached independently of the side rails for quick and easy access to both head and feet.



Trusted Solutions, Passionate People






Features and Options

Great comfort and stylish design

The *Invacare Octave* is designed specifically with the bariatric client in mind. The bed frame extends to the edges of the bed base for enhanced support and ease of transfer. An easy to use hand control allows the client to operate the bed's height and profiling function. The *Octave's* heavy duty motors enable a smooth and almost silent operation. The *Octave* is a substantial bed, however, its design includes all the well-known qualities of Invacare's Scandinavian design.

Delivery and storage

To make delivery as easy as possible, the *Octave* is delivered on a transportation rack. Dismounting and assembly of the *Octave* couldn't be simpler, despite the size of the bed. Furthermore, the transportation rack makes it possible to economise on storage space and makes the *Octave* very easy to move between locations.

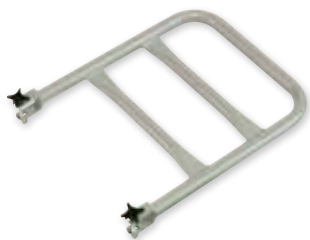

Support handle

For increased support whilst bed leaving.



Lifting pole

The handle can be adjusted in both height and depth, to offer support whilst re-positioning.



Extendable bedend

The bedend is extendable to offer mattress platform extension of 10 or 20 cm.

Technical data

<i>Invacare Octave</i>					
	Width outside: 119 cm	Length outside: 221 cm Length inside: 207 cm Mattress: 200 cm	Height adjustable: 40 - 80 cm Total height: 172 - 212 cm	Mattress support dimension: 81 - 34 - 33 - 52 cm (4-sectioned)	0° to 57° Backrest angle
<i>Invacare Octave</i>					
	0° to 30° Thigh angle	0 to 30° Legrest lift	Steel 180 kg Wood: 199 kg	SWL 385 kg/60 stone Max user weight: 350 kg / 56 stone	

3 year warranty

Complies with IEC 60601-2-52

Frame colour



Silver



Beech

NB:

The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Invacare Limited
Pencoed Technology Park -
Pencoed - Bridgend CF35 5AQ -
United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS -
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2016 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

Octave - UK - 03/2016 - 1487891



Yes, you can.®



For mattress compatibility, please refer to the user manual which can be found at www.invacare.co.uk