



*Invacare*<sup>®</sup> | **Propad<sup>™</sup> Mattress Overlay**

High specification pressure redistributing  
overlay



The **Propad Overlay** is designed to be placed on top of an existing conventional mattress, providing a high level of patient comfort with excellent pressure redistribution abilities. The overlay is lightweight, available in a range of sizes and suitable for those considered at 'High Risk'\* of developing pressure ulcers, making it ideal for the homecare environment.


## Effective Pressure Ulcer Prevention

The foam core has independent, geometrically-cut surface cells which help to reduce shear and friction, whilst increasing the overall support surface area, thereby reducing peak pressure.

## Even Weight Distribution

The **Propad Overlay** has larger, more supportive cells at the centre and smaller cells at each end, enabling the mattress overlay to achieve even weight distribution by conforming to the shape and form of the body and to cradle bony prominences.

## Features and Options

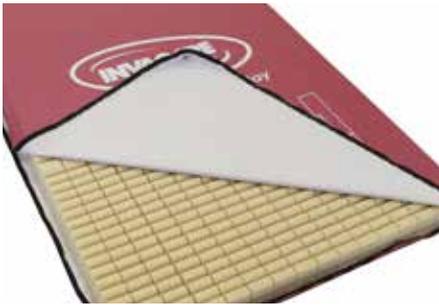

### Cover helps to alleviate shear and friction forces

Cover is multi-stretch, water-resistant and vapour-permeable. A white underside makes it easy to inspect to assist with infection control.



### Available in a range of sizes

Single - 187.5 x 87.6 x 8.5 cm,  
Double - 187.5 x 137.2 x 8.5 cm,  
King Size - 198 x 152.4 x 8.5 cm.



### Covers have an anti-slip base

Ensures safe use of the **Propad Overlay** and helps to prevent slippage, reducing the risk of falls.

## Technical data

For more information about this product, including the product's user manual, please visit [www.invacare.co.uk](http://www.invacare.co.uk).

	Size range	Height	Total product weight	Max. user weight	Max. cleaning temperature
<b>Propad Mattress Overlay</b>	Std sizes: 1875 x 876 mm <sup>1</sup> 1875 x 1372 mm <sup>1</sup> 1980 x 1524 mm <sup>1</sup> Additional sizes available on request	85 mm	5.5 kg	108 kg/ 17 Stone	80°

## Colour

Violet Quartz



1. Single
2. Double
3. King Size

## Warranty

3 Years

Manufactured to European Standards:  
BS 7177 Crib 5

\* Essential nursing care is pivotal in pressure ulcer prevention. This product will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value. The risk category in this instance is based on the Waterlow Scale.

Invacare Limited  
Pencoed Technology Park - Pencoed  
Bridgend CF35 5AQ - United Kingdom  
Tel: +44 1656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail:  
[ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS -  
County Dublin - Ireland  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [Ireland@invacare.com](mailto:Ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

Trusted Solutions,  
Passionate People



© 2017 Invacare International Sàrl. All rights reserved. All information quoted is believed to be correct at time of print.  
Invacare reserves the right to alter product specifications without prior consultation. Propad Overlay - UK - 10/2017 - 1504729