



Invacare[®]

Essential Plus

A cost effective pressure redistributing
mattress



The *Invacare Essential Plus* is an entry level pressure ulcer prevention mattress designed for the community healthcare environment. Practical and durable, this mattress offers comfort, support and pressure redistribution at a cost effective price. The *Essential Plus* is suitable for those considered at 'High Risk'¹ of developing pressure ulcers.

Effective Pressure Ulcer Prevention

The core of the *Invacare Essential Plus* is constructed from a series of 'keyhole cut' cells, which move independently to evenly distribute weight, promoting greater comfort and facilitating pressure ulcer prevention.

Cover is resistant to fluid ingress

The outer cover is a multi stretch, water resistant, polyurethane fabric providing excellent recovery and durability which helps to reduce shear and friction forces. The cover has a white inner liner allowing easy inspection for strikethrough or other cover damage.



Features and Options



Turn and rotate

Double sided castellated foam allows the mattress to be turned or rotated for increased product longevity.



Available in the following sizes:

L2000 x W855
Order code: 1576201
L1950 x W880
Order code: 1576200
L2050 x W880
Order code: 1629154








Zip around three sides

Allows easy inspection of the foam core and easy removal for laundering.


Technical data

For more information about this product, including the product's user manual, please visit www.invacare.co.uk.

Mattress sizes	Height	Total product weight	Max. user weight	Max. cleaning temperature	Warranty
 Essential Plus 2000 x 855 mm 1950 x 880 mm and 2050 x 880 x 150 mm	 140 mm	 14 kg	 140 kg	 80°	4 years foam 2 years cover Mattress is manufactured to: BS EN 597-1 Cig, BS EN 597-2 Match and BS 7177: Crib Source 5

Colour

Baltic Blue

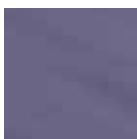

Warranty

4 years foam
2 years cover

- Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS -
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

Trusted Solutions,
Passionate People

