



Invacare®

COMET PRO

Power and performance in
a premium package



Whether you're popping to the local shops, taking longer trips into town, or looking to go further afield, the all new **Comet Pro** will take you the distance in comfort, safety and style.

Its powerful outdoor capabilities, combined with advanced safety features and options, ensures a speedy, powerful and reliable drive.

With its exceptional driving features and modern aesthetics, the **Comet Pro** offers power and performance in a premium package.



Premium performance and control

Thanks to its enhanced driving features and 2 x 75 Ah batteries, the **Comet Pro** allows you to enjoy longer and more challenging journeys beyond the city's limits. Its large 13" front and rear wheels, combined with the all new suspension system, allows the safe negotiation of kerbs and obstacles in a smooth and comfortable manner. Additionally, the automatic speed reduction will slow you down when driving around corners at high speed for additional stability and safety.

TÜV tested

The **Comet Pro** has been subject to independent German TÜV laboratory testing, giving you ultimate confidence and reassurance when driving. Advanced safety features, such as a two step disengaging lever and automatic speed reduction, offers additional peace of mind wherever you go.

Compact and modern design

Designed with both comfort and precision in mind, the **Comet Pro** offers exceptional manoeuvrability and control, even through tight and challenging spaces. Additionally, the new ergonomic steering ensures an optimal driving position, allowing for a more enjoyable and relaxing journey.



Features and options



Advanced suspension

All new suspension system for a smooth and comfortable drive over all terrains.



Easy drive direction stability

Keeps you in a straight line at high speeds for additional safety and stability.



Finite tiller movement

Easily adjustable tiller with stepless adjustments to enhance steering control. This coupled with the New head tube angle of 70 Degrees, help to give more responsive steering control.



Ergonomic steering

Our ergonomically designed steering system makes for easier, more responsive control.



Intuitive LCD Display

Easy to understand layout giving you complete confidence and control when driving.



Brake light

Clear and bright brake lighting warns people that the scooter is slowing down.




Disengaging lever

Two step disengaging lever prevents the scooter from accidentally free-wheeling.



Auto-speed reduction

Automatic speed reduction reduces speed around bends for additional safety.

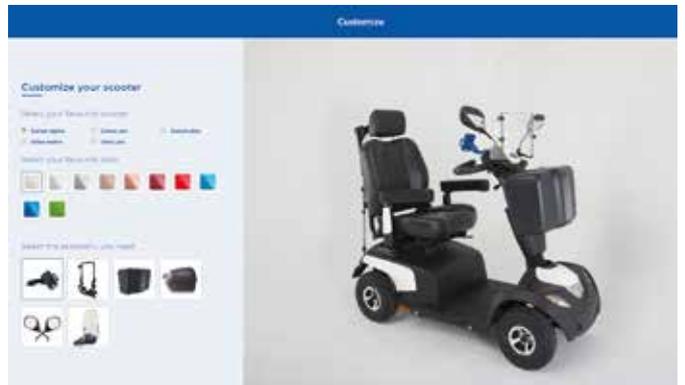


Customise your scooter with shrouds and accessories

Choose from 10 shroud colours and over 25 accessories for customising your scooter to your individual style. To find out more, visit our dedicated scooter website at: www.invacarescooter.com.

Technical data

For more comprehensive information about this product, including the product's user manual, please visit www.invacare.co.uk

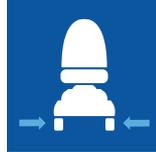	Backrest width	Seat depth	Seat to floor height	Backrest height	Driving unit width	Total length	Empty weight
<i>Comet Pro</i>	 510 mm	 460 mm	 Min. 440 mm Max. 510 mm	 630 mm	 665 mm - 685 mm	 1450 mm	 136 kg

	The heaviest part	Max. user weight	Battery capacity	Turning diameter	Max. climbable obstacle height	Ground clearance	Max. safe slope
<i>Comet Pro</i>	 approx. 63 kg	 160 kg (25st)	 2 x 12 V/75 Ah (C20) AGM	 2750 mm	 100 mm	 100 mm	 10°

	Max. driving distance	Max. driving speed	Motor capacity
<i>Comet Pro</i>	 55 km (34 miles)	 12.8 km/h (8MPH)	 550 W / 1500 W ¹

1. Peak performance

Colour of shrouds

Upholstery

Classic Colour Line

Marble White	Indian Blue	Jade Green	Diamond Silver	Jasper Red	Black Vinyl
					

Premium Colour Line

Carmine Red	Iconic Chrome	Electric Blue	Luxor Beige	Crystal Copper
				

Invacare Limited
 Pencoed Technology Park - Pencoed
 Bridgend CF35 5AQ - United Kingdom
 Tel: +44 1656 776 222
 Fax: +44 1656 776 220
 E-mail: uk@invacare.com
 Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
 Unit 5 - Seatown Business Campus
 Seatown Road - SWORDS -
 County Dublin - Ireland
 Tel. +353 1 8107084
 Fax +353 1 8107085
 Email: Ireland@invacare.com
www.invacare.ie

Trusted Solutions,
 Passionate People

