



Invacare[®] | **Octave**[®]

Setting new standards in Bariatric care



With a safe working load up to 60 Stone (385 kgs), the *Invacare Octave* is a robust yet stylish profiling bed solution that assists in the care and comfort of bariatric clients. With a fully electrical, 4-sectioned profiling mattress support as standard, the **Octave** provides optimal comfort and support to the client, as well as an ideal working environment for the carer. The frameless construction of the mattress platform allows improved side access to the client, and the bed ends can be detached independently of the side rails for quick and easy access to both head and feet.

Great comfort and stylish design

The **Octave** is designed specifically with the 'plus sized' bariatric client in mind. The bed frame extends to the edges of the bed base for enhanced support and ease of transfer. An easy to use hand control allows the client to operate the bed's height and profiling function. The **Octave's** heavy duty motors enable a smooth and almost silent operation. The **Octave** is a substantial bed, however, its design includes all the well-known qualities of Invacare's Scandinavian design.

Features and Options



Lifting pole

The handle can be adjusted in both height and depth, to offer support whilst re-positioning.



Extendable bedend

The bedend is extendable to offer mattress platform extension of 10 or 20 cm.



Support handle

For increased support whilst bed leaving.

Delivery and storage

To make delivery as easy as possible, the **Octave** is delivered on a transportation rack. Dismounting and assembly of the **Octave** couldn't be simpler, despite the size of the bed. Furthermore, the transportation rack makes it possible to economise on storage space and makes the **Octave** very easy to move between locations. Ideal for bed rental needs.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please visit www.invacare.co.uk.

	Width outside	Length outside Length inside	Height adjustment	Mattress support dimensions	Backrest angle	Thigh angle	Legrest lift	Total weight
<i>Invacare Octave</i>	1190 mm	2210 mm 2070 mm Mattress: 2000 mm	400 - 800 mm Total height: 1720 - 2120mm	810-340- 330-520 mm (4-sectioned)	0 - 57°	0 - 30°	0 - 30°	Steel 180 kg Wood: 199 kg


	Max. user weight
<i>Invacare Octave</i>	SWL 385 kg/ 60 stone 350 kg/ 56 stone



Complies with IEC 60601-2-52.

NB:
The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour

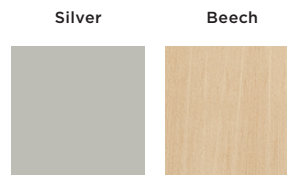

Warranty

3 years frame
3 years actuator
3 years control box
2 years hand control

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

