

Now with
fully welded
seams and a
waterproof zip

2 Years warranty
*High risk/
Very high risk

Invacare® | **Propad™ Renal Overlay**

High specification pressure redistributing overlay ideal for dialysis units



The *Invacare Propad Renal Overlay* has been developed in conjunction with clinicians specifically to provide additional comfort, effective pressure redistribution and security for patients undergoing dialysis treatment. Treatments, such as dialysis performed over a long period of time on a support surface that restrict free movement, can greatly increase the risk of skin damage as a result of direct pressure over bony prominences. Specifically cut castellations in zoned areas across the surface of the renal overlay not only enhance the pressure redistributing qualities of a dialysis chair, they will also address areas of increased vulnerability where there is a reduction or absence of subcutaneous fat (the heels or the elbows). Confidence for staff and patients is provided by fully adjustable securing straps to securely fit over most dialysis chairs, preventing movement of the overlay mattress during patient transfer. Infection control requirements can be easily met by the **recent addition of a fully sealed zip**, reducing the risk of fluid ingress to the internal components of the overlay. The vapour permeable wipe-able surface can be easily cleaned between patients to reduce valuable staff time and essential treatment delays.

Features and Options



Effective Pressure Ulcer Prevention

The high specification foam core has independent, geometrically-cut surface cells which help significantly reduce shear and friction. The **Renal Overlay** has larger, more supportive cells at the centre and smaller cells at the head and foot end. This enables improved weight distribution with the overlay able to conform to the shape of the patient's body.



Comfort for the Patient

Four securing straps, two of which are adjustable allow the overlay to profile when used in conjunction with an adjustable chair. The **Renal Overlay** is hinged at the foot, to allow perfect contouring to the vulnerable heel area with the benefit of increased comfort, when the foot plate is activated.



Size

213 x 57 x 8.5 cm



Specialist cover design

With infection control being paramount in busy renal units, the high frequency **welded seams and waterproof zip** with interlocking teeth on the renal cover are essential new additions to this overlay. The Renal overlay cover is multi stretch, water-resistant and vapour-permeable. This facilitates patient movement reducing the potential for tissue damage due to shear and friction.

Technical data

For more information about this product, including the product's user manual, please visit www.invacare.co.uk.

Colour

Warranty

Baltic Blue


2 years

Size range



Std sizes:
218 x 57 cm

Height

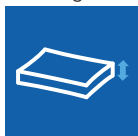

8.5 cm

Total product weight




5.5 kg

Max. user weight



108 kg

Max. cleaning temperature

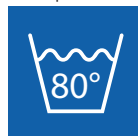

80°

Propad Mattress Overlay

Manufactured to European Standards: BS 7177 Crib 5

*Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgment and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgment). The risk category in this instance is based on the Waterlow Scale

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS -
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

INVACARE
Yes, you can.®



© 2017 Invacare International Sàrl. All rights reserved. All information quoted is believed to be correct at time of print.
Invacare reserves the right to alter product specifications without prior consultation. Propad Renal Overlay - UK - 10/2017 - 1585805