



Invacare®

Reliant 350 Stand Assist

For safe and comfortable, single carer assisted transfers



Features and Options

The *Invacare Reliant 350*

Stand Assist is a compact and maneuverable mobile lifter that provides safe and comfortable assisted transfers for those with limited mobility or rehabilitation needs. The *Invacare Reliant 350* can be used to transfer a client that cannot transfer unaided, but has some level of weight bearing ability, from one seated position to another. The robust design with a weight limit of 158 kg allows for additional support whilst transferring, preventing unnecessary strain on the carer, making it a suitable option for transfers carried out by a single carer*. • 5 year warranty: Frame, Actuator, Control box

• 2 year warranty: Hand control




Independent Charger

An off-board charger and spare battery allow for 24 hour use. Please note: The battery cannot be charged on the hoist.


Stand Assist Sling

Designed to be used with a **Stand Assist** lifter, this sling features the latest innovative fabric for improved comfort and handling. Available in sizes; S, M, L & XL.

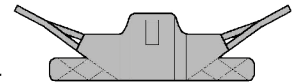

Transfer Stand Assist Sling

Designed for **Stand Assist** lifters, for those who require additional leg support, this sling provides a supportive and comfortable transfer. Available in sizes; S, M, L & XL.

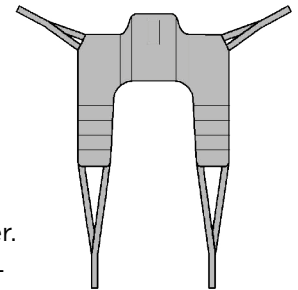

Adaptable and easy to use

Reliant 350 features a multi-adjustable lift boom and a selection of sling mounting points to provide an adjustable and ergonomic lift angle. Large rear and small front castors help manoeuvre the **Reliant 350** into small or confined spaces and the base has a lever adjustment mechanism for added stability. In addition, it has a generous non slip footplate providing a stable surface for the feet, as well as a height-adjustable padded knee guard.

Safe and reliable

The **Reliant 350** provides a range of safety features including an emergency stop button, emergency mechanical descent system, as well as an anti-entrapment device which stops the boom from lowering if it meets resistance. The battery recharge indicator is visual and audible. For transport and storage, the **Reliant 350** can be easily dismantled.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local **Invacare** website.

Max. height	Min. height	Overall length of base	Inside width: legs opened	Outside width: legs closed	Height of legs	Number of cycles	Total product weight
1670 mm	1020 mm	900 mm	940 mm	660 mm	110 mm	100 - 200 lifts per full charge	49 kg

Invacare
Reliant 350

Max. user weight



158 kg

Invacare
Reliant 350

Frame colour

Beige



Warranty

Frame, control box & actuator	5 years
Hand control	2 years
Battery	1 year

The *Invacare Reliant 350 Stand Assist* is part of the 'Safe Patient Handling' range.



safepatienthandling

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

