



*Invacare*<sup>®</sup> | **Softform<sup>®</sup> Premier Visco**

Pressure redistributing visco elastic foam  
mattress



The *Invacare Softform Premier Visco* is a high specification pressure redistributing mattress suitable for those considered at 'High Risk'<sup>1</sup> of developing pressure ulcers. The mattress features a high density visco elastic (memory) foam insert, which is temperature-sensitive to the body, allowing it to conform to body shape. This provides exceptional comfort and facilitates pressure redistribution by minimising peak pressures. The '**Strikethrough Resistant Technology**'<sup>2</sup> cover has been designed to manage the clinical impact in relation to pressure care and infection control with regards to risk of cross infection and financial impact, including time and resource management.


### Comfort and pressure redistribution combined

The *Invacare Softform Premier Visco* is constructed from a range of superior quality materials. The lower 'U' shaped layer is made up of robust high density foam, to prevent the mattress bottoming out and improve longevity. The central core includes a high density, visco elastic (memory) foam top section which warms to conform to body shape, maximising the body contact area to give good pressure reduction.

### Strikethrough Resistant Technology® (SRT)

The cover is made from a highly chemically resistant polyurethane polymer making it much more resistant to abrasion and damage and swells less with water contact. When cleaning a cover with '**Strikethrough Resistant Technology**'<sup>®</sup>, the material properties changes less than standard fabrics and thus reverts to its natural state very quickly. The time in which the cover is more prone to physical damage is significantly reduced.

## Features and Options

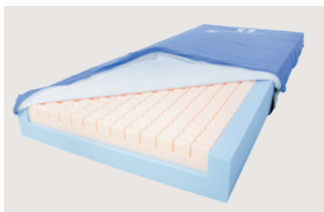


#### Does not require turning or rotating

This mattress does not require turning or rotating, significantly minimising manual handling issues.



#### Modular construction

Enables easy replacement of parts, increasing product longevity. Insert pad can be interchanged with standard castellated foam, if required.

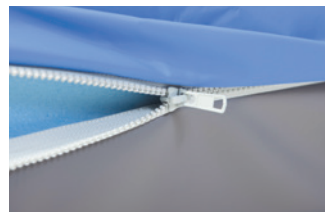

#### Full-length flap conceals zip

The full length flap protects the zip. This reduces the possible risk of fluid ingress, preventing contamination of the foam core.



#### Added protection

Multi-stretch, breathable, waterproof cover with corner welds, reducing the opportunity for accumulation of dust and liquid residues.

## Technical data

For more information about this product, including the product's user manual, please visit [www.invacare.co.uk](http://www.invacare.co.uk).

### Colour

	Width	Length	Height	Total product weight	Max. user weight	Max. cleaning temperature	Colour
<b>Softform Premier Visco</b>	660 - 900 mm*** 910 - 1370 mm** 1380 - 1525 mm*	1410 - 2080 mm	152 mm	13 kg	247.6 kg/ 39 Stone	80°	Violet SRT
	*King size, **Double size, ***Single size	Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match, BS7177: 1996 Crib Source 5		<b>Warranty</b> Cover 4 years Foam 4 years			



1. Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

2. Please note: The cover is not Strikethrough proof. Misuse or incorrect care and handling can cause damage to the mattress cover.

The Double and Kingsize mattresses are not available with the Strikethrough Resistant Technology<sup>®</sup> cover.

**Invacare Limited**  
Pencoed Technology Park - Pencoed  
Bridgend CF35 5AQ - United Kingdom  
Tel: +44 1656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail:  
[ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

**Invacare Ireland Ltd**  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS -  
County Dublin - Ireland  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [Ireland@invacare.com](mailto:Ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

