

matrix[®]
SEATING SERIES

Flo-tech™ Solution



Slimline cushion with unique dual
fluid sac system



The *Invacare **Matrix Flo-tech Solution*** cushion uses a combination of a slim, shaped foam cushion base and a dual layer of fluid sacs over the ischial and sacral area. The dual fluid sacs and the unique moulded foam shape of the solution increase the support surface area. This delivers a superior level of protection for patients identified as being at Very High Risk of developing pressure ulcers, or for those with already established tissue damage.

Features and Options



Ingress Resistant Cover

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle.

Available in 44 sizes

The *Invacare **Matrix Flo-tech Solution*** is available in a range of 44 sizes, from 38 x 41 to 51 x 46 cm (38 x 41 to 510 x 460 mm). (See size chart for details).

Invacare *Matrix Flo-tech Deep Solution* option

Option of a 10 cm (100 mm) "Deep" foam base is also available in 12 sizes. (See size chart for details).

Ultra-fresh™ added to all cushions for antimicrobial and odour protection

Antimicrobial treatment that has been added to all ***Flo-tech*** foam to help make it resistant to the growth of odour causing bacteria. Helps provide longer lasting freshness.

Slimline contoured foam base with large dual fluid sac covering the entire top surface of the cushion







The contact area between the user and the foam cushion is protected, as the fluid displaces to redistribute pressure away from bony prominences. The additional layer of fluid is incorporated under the ischial and sacral area to further enhance the pressure reducing qualities of this cushion, making this cushion suitable for Very High Risk users.

Slimline contoured cushion

Some users prefer a slimline cushion, finding it more compatible with their daily living activities. The low profile is particularly beneficial for the facilitation of side transfers.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.


| | Various sizes | Height | Total product weight | Max. user weight ¹ | Max. cleaning temperature |
|---|---|---|---|--|--|
|  |  |  |  |  |  |
| Matrix Flo-tech Solution | A range of 44 sizes | 78 mm | 3.08 kg | 140 kg | 80° |

Manufactured to European Standards: BS EN 1021-1 Cig, BS EN 1021-2 Match and BS 7176 (Crib 5)

Disclaimer

1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.

Colours

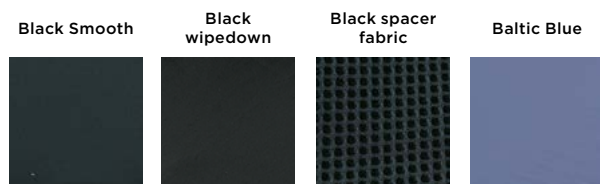

QR Code

Invacare Ltd
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel +44 1 656 776 222
E-mail: uk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS, County Dublin - Ireland
Tel. +353 1 8107084
Email: Ireland@invacare.com
www.invacare.ie

