

# ***Reliant 350 Stand Assist***



For safe and comfortable,  
single carer assisted transfers



The *Invacare Reliant 350 Stand Assist* is a compact and maneuverable mobile lifter that provides safe and comfortable assisted transfers for those with limited mobility or rehabilitation needs. The *Invacare Reliant 350* can be used to transfer a client that cannot transfer unaided, but has some level of weight bearing ability, from one seated position to another. The robust design with a weight limit of 158 kg allows for additional support whilst transferring, preventing unnecessary strain on the carer, making it a suitable option for transfers carried out by a single carer\*. • 5 year warranty: Frame, Actuator, Control box

- 2 year warranty: Hand control

### Adaptable and easy to use

*Reliant 350* features a multi-adjustable lift boom and a selection of sling mounting points to provide an adjustable and ergonomic lift angle. Large rear and small front castors help manoeuvre the *Reliant 350* into small or confined spaces and the base has a lever adjustment mechanism for added stability. In addition, it has a generous non slip footplate providing a stable surface for the feet, as well as a height-adjustable padded knee guard.

## Features and Options



### Independent Charger

An off-board charger and spare battery allow for 24 hour use. Please note: The battery cannot be charged on the hoist.


### Stand Assist Sling

Designed to be used with a *Stand Assist* lifter, this sling features the latest innovative fabric for improved comfort and handling. Available in sizes; S, M, L & XL.



### Transfer Stand Assist Sling

Designed for *Stand Assist* lifters, for those who require additional leg support, this sling provides a supportive and comfortable transfer. Available in sizes; S, M, L & XL.

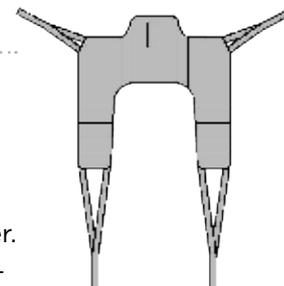

### Safe and reliable

The *Reliant 350* provides a range of safety features including an emergency stop button, emergency mechanical descent system, as well as an anti-entrapment device which stops the boom from lowering if it meets resistance. The battery recharge indicator is visual and audible. For transport and storage, the *Reliant 350* can be easily dismantled.

## Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please visit [www.invacare.co.uk](http://www.invacare.co.uk).

| Max. height | Min. height | Overall length of base | Inside width: legs opened | Outside width: legs closed | Height of legs | Number of cycles                | Total product weight |
|-------------|-------------|------------------------|---------------------------|----------------------------|----------------|---------------------------------|----------------------|
| 1670 mm     | 1020 mm     | 900 mm                 | 940 mm                    | 660 mm                     | 110 mm         | 100 - 200 lifts per full charge | 49 kg                |

Invacare  
**Reliant 350**

| Max. user weight |
|------------------|
| 158 kg           |

Invacare  
**Reliant 350**

### Frame colour

Beige



### Warranty

|                               |         |
|-------------------------------|---------|
| Frame, control box & actuator | 5 years |
| Hand control                  | 2 years |
| Battery                       | 1 year  |

The *Invacare Reliant 350 Stand Assist* is part of the 'Safe Patient Handling' range.



safe patient handling

Invacare Limited  
Pencoed Technology Park - Pencoed  
Bridgend CF35 5AQ - United Kingdom  
Tel: +44 1656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail: [ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS - County Dublin - Ireland  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [Ireland@invacare.com](mailto:Ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

