



Compatible with the
NEW Active 2 RX pump

Invacare® | **Softform® Premier Active Cushion**

A comfortable cushion that supports 24hr
pressure care management



The *Invacare Softform Premier Active Cushion* has been designed for patients who are considered at 'High Risk' of developing pressure ulcers and require pressure care management 24 hours a day. The unique design features an alternating air insert beneath the castellated foam of the *Softform Premier Active Cushion*, transforming the cushion to a dynamic surface in seconds. This is designed to operate from the same pump as the *Softform Premier Active 2* Mattress ensuring the patient has continuity of care with 24 hour pressure care management.

An innovative active cushion

The air cells work on a two up, two down function beneath the insert, delivering a very low impact at the point of change, thereby delivering high levels of comfort in addition to improved pressure redistribution.

How it works:



Unplug *Softform Premier Active 2* mattress hose from *Softform Active 2 RX* pump.



Plug hose from *Softform Premier Active* cushion into *Softform Active 2 RX* pump.



Switch on pump.



Patient can now sit on an alternating surface, to provide 24hr protection.

Features and Benefits



Compatible with *Softform Premier Active 2* Mattress

The same pump can be used for both mattress and cushion, offering a 24 hr support system.






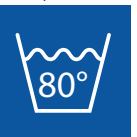


All new *Active 2 RX* pump

Pressure can be maintained at a level appropriate to individual patients, offering continuity of care.

Technical data

For more information about this product, including the product's user manual, please visit www.invacare.co.uk.

Width	Length	Height	Total product weight	Max. user weight	Max. cleaning temperature	Warranty
						Cover: 12 months Foam: 12 months Pump: 24 months Air Insert: 12 months
<i>Softform Premier Active Cushion</i>	460 mm	510 mm	152 mm	1.5 kg	17 stone (108 kg)	80°

Colour*

Baltic Blue



Product codes

Standard Active Cushion	1523384
Air Insert	1523383
Active 2 RX Pump	1627629

For information on the *Softform Active 2 RX* pump, please see the *Softform Premier Active 2* brochure. Part number: 1552287

Manufactured to European Standards:
Foam: BS 5852: Crib 5; Cover BS7175: Crib 5.

*Please note, the colour of the cushion may vary to that shown in this brochure

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS -
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

