


# Softform<sup>®</sup> Premier Visco

Pressure redistributing visco elastic mattress

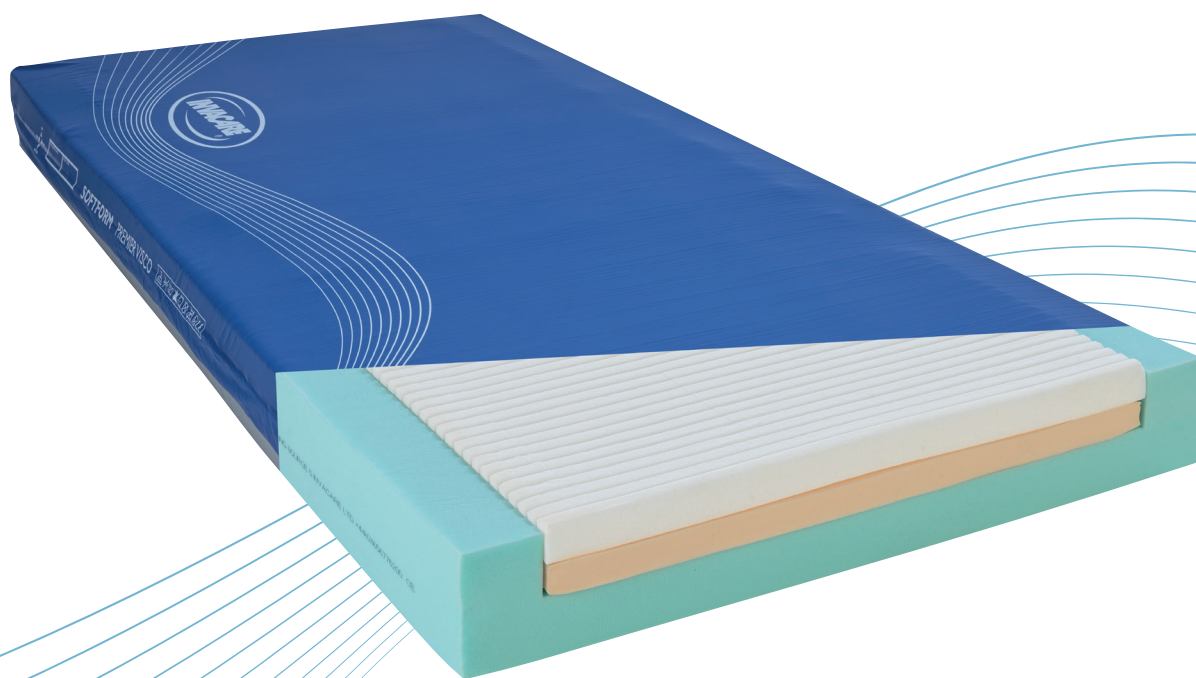

Yes, you can.<sup>®</sup>

# Effective pressure injury prevention

The Invacare Softform Premier Visco is a pressure redistributing mattress featuring a deep visco insert pad, which is high-density and temperature-sensitive to the body, allowing it to conform to the user's shape. This affords exceptional comfort and facilitates pressure reduction by minimising peak pressures.

The temperature-sensitive visco elastic insert pad within the Invacare Softform Premier Visco warms to conform to the user's shape, maximising the body


contact area and thereby facilitating good pressure reduction. The cover is a two-way stretch polyurethane material with high frequency welded seams; and a full length flap conceals the zip.

The high quality visco elastic foam featured in this mattress delivers an exceptionally comfortable support surface.



## Cover is resistant to fluid ingress

- › Features include two-way stretch **waterproof polyurethane cover** with high frequency welded seams and a full length flap concealing the zip.
- › Also has a toughened PU coated base to prolong the longevity of the mattress.

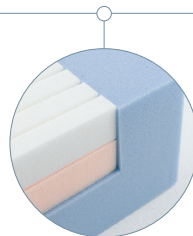

## Modular Construction

- › Enables easy **parts replacement**.
- › Helps to increase product lifetime.



## Orientation symbol

- › The orientation symbol ensures that the mattress is **positioned correctly** in the bed.

## Technical data ▼

For more information about this product, including the user manual, please visit our website: [www.invacare.co.uk](http://www.invacare.co.uk)



Width ▼



Height ▼



Length ▼



Weight ▼



User weight ▼



Clean temperature ▼



Baltic Blue

Colour

Single size:  
660 - 900 mm  
Double size:  
910 - 1370 mm  
King size:  
1380 - 1525 mm

152 mm

1410 - 2080 mm

13 kg

247.6 kg

95 °C

Manufactured to European Standards:  
BS EN 597-1 Cig,  
BS EN 597-2 Match.  
Foam Interior:  
BS 5852: 1990 Crib Source 5; Cover BS7175:  
1989 Crib Source 5.

1. Essential nursing care is pivotal in pressure injury prevention. This mattress will positively contribute to the outcome of a pressure injury prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure injuries. A range of assessment scales can be used as a formal method of assessing risk from pressure injury development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

T: +44 1656 776222  
[uk@invacare.com](mailto:uk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)  
[www.invacare.ie](http://www.invacare.ie)

Invacare UK & Ireland  
Unit 4 Pencoed Technology Park  
Bridgend CF35 5AQ



Yes, you can.

