


Softform[®] Excel

Pressure redistributing mattress

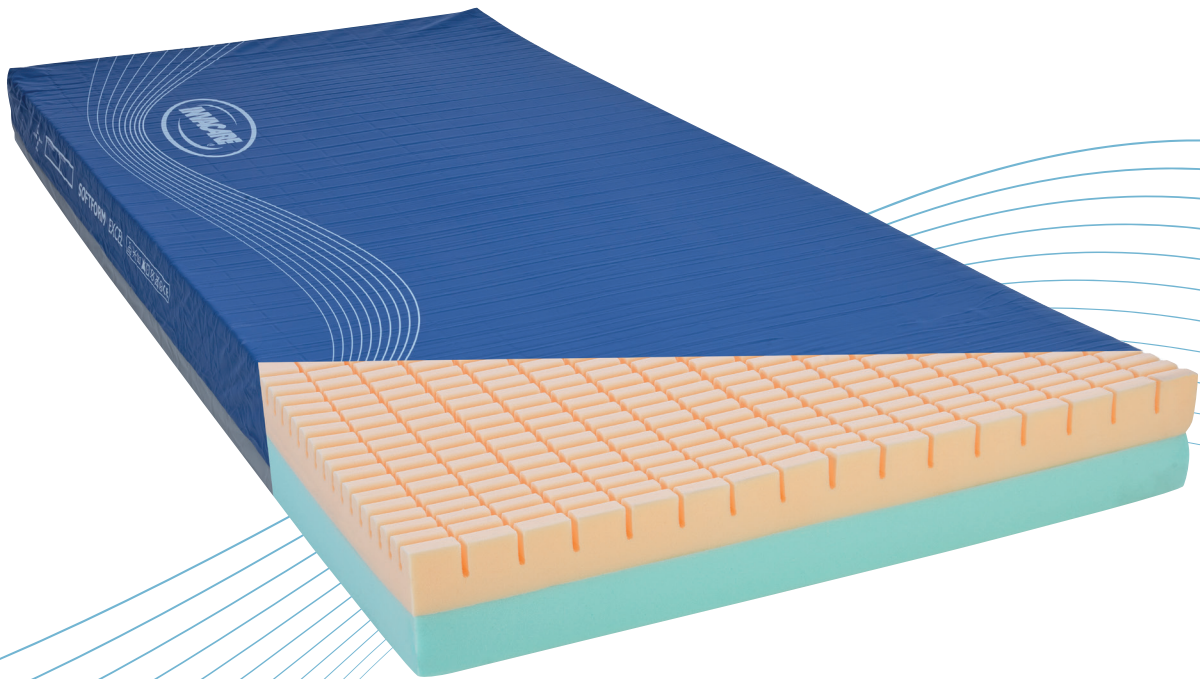

Yes, you can.[®]


Effective pressure injury prevention

The Invacare Softform Excel provides a high level of comfort and is constructed from high specification foam which provides even weight distribution for greater comfort and pressure redistribution. With a robust foam base, it prevents the mattress from bottoming out and is suitable for those considered at risk of developing pressure injuries.

The core of the Invacare Softform Excel is constructed from a series of 'keyhole cut' cells, which move independently to evenly distribute weight, promoting greater comfort and facilitating pressure redistribution.

The robust base features a firm layer of high density foam to provide substantial support.

The cover fabric provides excellent recovery and durability, and helps to reduce shear and friction forces. The base of the mattress is a tough coated polyurethane (PU) which protects the mattress when used with electric profiling bed frames.

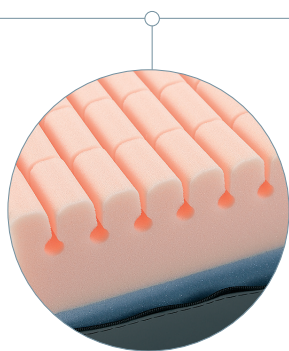

Strong, durable construction

- › The Invacare Softform Excel is an **established choice** of healthcare professionals.



Zip around three sides

- › Allows **easy inspection** of the foam core.



Orientation symbol

- › The orientation symbol ensures that the mattress is **positioned correctly** in the bed.

Technical data ▼

For more information about this product, including the user manual, please visit our website: www.invacare.co.uk



Width ▼



Height ▼



Length ▼




Weight ▼



User weight ▼



Clean temperature ▼

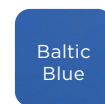

Baltic Blue

Colour

Single size:
800, 855, 880,
900 mm

152 mm

1950 mm

10 kg

140 kg

95 °C

Manufactured to European Standards:
BS EN 597-1 Cig,
BS EN 597-2 Match.
Foam Interior:
BS 5852: 1990 Crib Source 5; Cover BS7175:
1989 Crib Source 5.

1. Essential nursing care is pivotal in pressure injury prevention. This mattress will positively contribute to the outcome of a pressure injury prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure injuries. A range of assessment scales can be used as a formal method of assessing risk from pressure injury development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

Invacare Ltd
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1656 776 222
E-mail: uk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus Seatown Road -
SWORDS - County Dublin - Ireland Tel. +353 1
8107084
Email: Ireland@invacare.com
www.invacare.ie



Yes, you can.