

Softform Odstock Wedge



Pressure reducing wedge



The *Invacare Softform Odstock Wedge* was developed in partnership with clinicians to assist in the reduction of interface pressures on the sacrum and heels without significantly increasing pressure on the thighs or calves. It is designed to improve and maintain patient posture in the recumbent position, and to reduce the tendency for the user to slip down the bed, thereby reducing shearing forces on the sacrum. The **Odstock Wedge** is ideal for improving posture in those with back discomfort.



Features and Options




High frequency welded seams

Ensures core is protected from fluid ingress, increasing the product integrity.



High quality construction

The two way stretch cover, is waterproof and vapour permeable. The foam core is constructed from, high density foam.

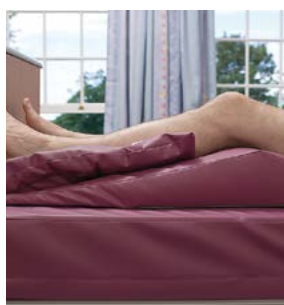

Designed to be used in conjunction with the *Invacare Heelpad*

Used in conjunction, these two products combine to significantly reduce pressure on the vulnerable heel area.

Reduces pressure on the heels







Clinical research has effectively demonstrated that the degree of knee flexion necessary to achieve substantial levels of pressure relief is relatively small. The *Invacare Softform Odstock Wedge* was designed by clinicians to reduce pressure on the heels through flexing the knee and hip and elevating the knee.

Reduces lumbar lordosis

By tilting the pelvis, the **Odstock Wedge** also reduces lumbar lordosis – desirable following back surgery and ideal for improving the sleeping posture of patients with back problems.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local **Invacare website**.

Width	Length	Height	Total product weight	Max. cleaning temperature	Warranty	
						Manufactured to European Standards: Interior Foam: BS 5852: Crib 5, Cover: BS7175: Crib 5
570 mm	865 mm	115 mm	1.5 kg	80°	1 year	

Refer to instructions for use and relevant infection control guidelines.

Colours

Blue



T: +44 1656 776222
uk@invacare.com
www.invacare.co.uk
www.invacare.ie

Invacare UK & Ireland
Unit 4 Pencoed Technology
Park
Bridgend CF35 5AQ

