

AQUATEC®

Ocean Ergo

adjustments and accessories



Simple adaptations that make a difference

Adjustable, ergonomic and versatile



Yes, you can.®

Adjustments and accessories

Our Ocean Ergo range offers an advanced solution in shower chairs. Each chair is designed to offer high levels of adjustments and configurability as standard. For further customisation, our latest range of accessories provides solutions for extra comfort, support and safety.



Footplates

- > Angle adjustable, flip-up and textured footplates with heel loops
- > No tools required

Footrests

- > Height adjustable, swing away and detachable footrests
- > No tools required

Armrests

- > Ergonomic, swing away armrests can be adjusted to different heights

Seat height

- > Seat to ground height can be adjusted between 440-600 mm* or 475-600 mm**

*Ocean Ergo, Ergo XL and 24" Ergo models
**Ocean 24" Ergo XL and VIP models

- > No tools required

Accessories

Seating and backrests



Soft seat

Redesigned to mirror the contours of the new ergonomic seat plate, the softer feel gives users with sensory or skin protection issues added relief. Smaller aperture also available for shorter or thinner users. Easy to fit and clean.



Soft seat insert

Available as a standard version or smaller version. Can be used with the soft seat when the commode function is not required.



Ergonomic soft seat

Manufactured from viscoelastic foam, it offers a much softer feel for those with skin integrity or sensory issues. The foam creates greater depth of immersion promoting increased pelvic stability and positioning.



Variable soft seat

Hygiene cut out can be positioned in all four directions, improving access for carers assisting with personal hygiene.



Universal soft seat

Comfortable, closed seat surface for when commode function is not in use.



Solid backrest

The solid back support offers an alternative to the standard tension adjustable backrest, when firmer support is required to achieve a more upright posture. Easy to clean.



Armrest pads

New design offers improved support and comfort for the user (not for XL armrests).



XL back

Increases the width between the armrests by 80 mm. Can be retrofitted to standard Ocean Ergo models.



Wider XL armrests

Increases the distance between the armrests by 40 mm on each side.



Soft backrest cushion

For increased comfort. Can also be used to reduce the overall seat width and depth, helping to position smaller users.



Comfort backrest

An alternative to the standard backrest with a softer feel material for increased comfort.

Safety



Safety bar

For additional safety and reassurance.



Safety lock for armrest

Locks armrest into the vertical or horizontal position.



Front and rear anti-tippers

For additional safety plus front and rear stability.



Tilt-stop

An aid for those who are hoisted in/out of the chair, the tilt-stop screws ensure the frame remains horizontal.

Positioning



Safety belts

Broad pad manufactured from soft material with single pull buckle. Can be used across the pelvis or in addition across the chest for positioning and safety. Easy to attach. Machine washable.




Lateral supports

Adjustable in height via three possible positions on the backrest cane. Broad lateral pads can be repositioned for depth adjustment, comfort and effective pressure redistribution.



Ergonomic headrest holder

Offers multidirectional adjustment to correct or accommodate for more complex head/neck positioning needs. **Pad not included.**

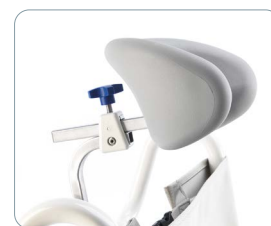

Ergonomic headrest/neckrest

Offers exceptional management of the head position. **Neckrest shown.**

Toileting



Calf and amputee support with foot support

Offers both height and angle adjustment with a contoured pad for added positioning and comfort. Including a foot support for those who may need to elevate their lower limb(s) for medical reasons or comfort.



Calf and amputee support

Excluding footrest.



Splash guard

Clips on easily for use with the standard plastic seat, standard soft seat and small soft seat.



Toilet pan

Ergonomic handle allows the pan to be easily removed and replaced. Features an integrated lid holder, so the lid is not misplaced during use.

Customise each of our Ocean Ergo chairs

Our Ocean Ergo range of shower chairs can all be customised to suit a wide number of needs. By making a simple adjustment or by adding one of our easy-to-fit accessories, each of our modular shower chairs are more versatile than ever before.



Ocean Ergo

An entry level attendant propelled model



Ocean Ergo XL

An entry level attendant propelled model with a heavier maximum user weight



Ocean 24" Ergo

An entry level self-propelling model



Ocean 24" Ergo XL

An entry level self-propelling model with a heavier maximum user weight



Ocean VIP Ergo

A tilt-in-space model for additional seating and positioning support



Ocean Dual VIP Ergo

A tilt-in-space model with recline for additional seating and positioning support

Compatibility

✓ Compatible
x Not compatible
+ Included

Most accessories are cross compatible with all Ocean Ergo models.

	Product code	Ocean Ergo	Ocean 24" Ergo	Ocean Ergo XL	Ocean 24" Ergo XL	Ocean VIP Ergo	Ocean Dual VIP Ergo
Standard soft seat	AP1603142	✓	✓	✓	✓	✓	✓
Soft seat insert	AP1603143	✓	✓	✓	✓	✓	✓
Small soft seat	AP1603145	✓	✓	✓	✓	✓	✓
Small soft seat insert	AP1603152	✓	✓	✓	✓	✓	✓
Splash guard	AP1603155	✓	✓	✓	✓	✓	✓
Universal soft seat	AP1603154	✓	✓	✓	✓	✓	✓
Ergonomic special soft seat	AP60121311	✓	✓	✓	✓	✓	✓
Variable soft seat	1544469	✓	✓	✓	✓	✓	✓
Safety bar 2.0	AP60123732	✓	✓	x	x	✓	✓
Safety bar 2.0 XL (only with armrests 2.0)	AP60123771	x	x	✓	✓	x	x
Safety bar XL (only for XL chairs with old armrests)	1470078	x	x	✓	✓	x	x
Chest belt	1470081	✓	✓	✓	✓	✓	✓
Pelvic belt	1470082	✓	✓	✓	✓	✓	✓
XL back	AP1604015	✓	✓	+	+	x	x
Solid backrest	1575299	✓	✓	x	x	✓	✓
Soft backrest cushion	1558769	✓	✓	✓	✓	✓	✓
Armrest pad (1 pad; only for new armrests 2.0)	AP1643748	✓	✓	✓	✓	✓	✓
Wider XL armrests	1470086	✓	✓	✓	✓	✓	✓
Safety lock for XL armrest	1470077	✓	✓	✓	✓	✓	✓
Safety lock 2.0 (for standard armrests 2.0)	AP1643760	✓	✓	✓	✓	✓	✓
24" wheels (set)	16379	✓	+	✓	+	x	x
Calf and amputee support right	AP1643573	✓	✓	✓	✓	✓	✓
Calf and amputee support left	AP1643574	✓	✓	✓	✓	✓	✓
Calf support with foot support right	AP1643543	✓	✓	✓	✓	✓	✓
Calf support with foot support left	AP1643572	✓	✓	✓	✓	✓	✓
Lateral support	1535077	✓	✓	✓	✓	✓	✓
Ergonomic headrest holder	AP1531613	✓	✓	✓	✓	✓	✓
Ergonomic headrest	1531465	✓	✓	✓	✓	✓	✓
Ergonomic neckrest	1531608	✓	✓	✓	✓	✓	✓
Toilet pan	10230	✓	✓	✓	✓	✓	✓
Toilet pan holder	AP1602078	✓	✓	✓	✓	✓	✓
Front and back anti-tipper Ergo VIP models	AP1603992	x	x	x	x	✓	✓
Front anti-tipper Ocean Ergo models	AP1601986	✓	✓	✓	✓	x	x
Comfort soft backrest VIP Ergo models	AP1603606	x	x	x	x	✓	✓
Comfort soft backrest Ocean Ergo models	AP1604136	✓	✓	x	x	x	x
Comfort soft backrest Ocean Ergo XL models	AP1604016	x	x	✓	✓	x	x
Tilt-stop screw set	AP60129156	x	x	x	x	✓	✓

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie