



**Yes, you can.<sup>®</sup>**

## **Invacare<sup>®</sup> Platinum<sup>®</sup> 9**

### **Reduce costs with the *Invacare<sup>®</sup>* high flow oxygen concentrator**

The *Invacare<sup>®</sup> Platinum<sup>®</sup>* range of oxygen concentrators offers the best value added and reliable solution to the treatment of clients requiring long term oxygen therapy.

The latest addition to the range is the **Platinum<sup>®</sup> 9** high flow concentrator which is fully compatible with the *Invacare<sup>®</sup> HomeFill<sup>®</sup>* system.

Whilst optimising oxygen concentration, the **Platinum<sup>®</sup> 9** substantially reduces the costs associated with maintenance and delivery. This concentrator enables clients to live at home in comfort and safety where they can appreciate it's simplicity of operation including features such as the integrated alarms which alert any changes in equipment performance.



## ▶ High performance, reliability and simplicity

Powerful and reliable, the *Invacare® Platinum® 9* is very user friendly. With a 9 litre oxygen capacity and 95,6 – 87% purity the *Platinum® 9* is suitable for clients requiring high levels of oxygen supply. The unit is easy to use and manoeuvre with clear concise controls, integrated handle and compact castors. Other features include:

- Easy access to the flow indicator
- On board O<sub>2</sub> sensor monitors the level of oxygen
- On and off switch
- Automatic circuit breaker
- Easy access to the humidifier bottle
- Twin head compressor ensuring better reliability



## ▶ Value added associated with cost control

The *Platinum® 9* has been designed to have the minimal number of parts so that low level of critical spares need to be retained - very cost effective! With no delivery, low maintenance costs together with a five year warranty, the *Platinum® 9* gives you the best value added concentrator.



**Secured power cord connection**



**Only one filter to change when dirty**




**HomeFill® connection port**

## Technical specifications

### Performance

Adjustable flow	1 – 9 l/min	Oxygen concentration	
Oxygen outlet pressure	62 kPa	2,7 l/min	95,6 – 93%
Oxygen outlet temperature	+2,2 / ambient temperature	8 l/min	95,6 – 91%
Sound level	< 50 db(a)	9 l/min	95,6 – 87%
Battery	no		
Low oxygen flow sound alarm	< 1 l/min		
Concentration alarms	> 85% Green light		
	> 73% Yellow light		
	< 73% Red light and sound		
Self-diagnostic system	failure code		


### Technical data

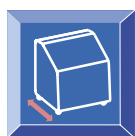

36,5 cm



46,7 cm

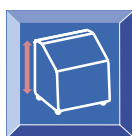

67 cm



24 kg



230 V, 50-60 Hz



475 W

All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation.



INVACARE® Limited  
 Pencoed Technology Park · Pencoed · Bridgend CF35 5HZ  
 Tel: +44 1 656 776 222 · Fax: +44 1656 776 220  
 E-mail: uk@invacare.com · [www.invacare.co.uk](http://www.invacare.co.uk)  
 Sales Order E-mail: ordersuk@invacare.com

INVACARE® Ireland Ltd  
 Unit 5 · Seatown Business Campus  
 Seatown Road · SWORDS · County Dublin  
 Tel. +353 1 8107084 · Fax +353 1 8107085  
 Email: Ireland@invacare.com · [www.invacare.ie](http://www.invacare.ie)

