


Invacare® *Flo-Shape cushion*



The *Invacare Flo-shape cushion* is a **103 mm** high cushion suitable for those **'At Risk'** of developing pressure ulcers.

Used predominantly on the Invacare Rea® range of wheelchairs, the cushion allows for **pressure redistribution** and **reduces peak pressure** to the most vulnerable areas.

The surface of the *Flo-shape cushion* is **contoured** into sections through a **mould** construction, giving a **three-dimensional** form that cannot be emulated in cut foam. This allows for weight to be **dispersed** away from key pressure points.

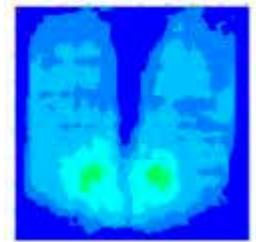


The anatomical shape of the cushion:

- **naturally supports** the ischial tuberosity, reducing the tendency for the person to slide forwards in the seat which can cause shear and friction
- provides a **stable surface** that **aids central pelvic placement** and **improved thigh alignment**, making it suitable for those with reduced sitting balance, such as a neurological condition.

Quick facts about the cushion

- Moulded foam not cut foam
- Contoured top surface
- Comfortable, stable and redistributes pressure
- Suitable for bony or heavier people
- Two-way stretch, water resistant and vapour-permeable Dartex cover
- Black smooth (plush) and black spacer fabric covers also available
- Available in a range of sizes
- Maximum user weight of 190.5 kg

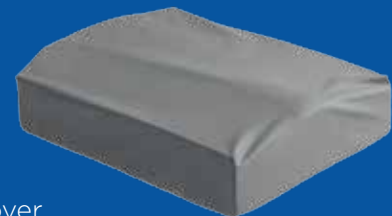

Technical Details

Foam	Grade Ref & Colour: FG 70/150 Grey Nominal Density (Kg/M): 65 – 75 Nominal Hardness (N): 150 Fatigue Class: Very Severe
------	--

Fire Testing	EN 597 - 1 & 2 BS 7176: 1995 crib 5
--------------	--

* Maximum user weight limit depends on cushion size

Knowledge bank

Invacare Limited
Pencoed Technology Park -
Pencoed - Bridgend CF35 5AQ
- United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business
Campus
Seatown Road - SWORDS -
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

rea
Confidence in posture solutions

INVACARE
Yes, you can.®