

matrix[®]
SEATING SERIES

Flo-tech™ Image



Slimline cushion with fluid sac



The *Invacare* **Matrix Flo-tech Image** cushion uses a combination of a slim, shaped foam cushion base and a fluid sac covering the sacral area, increasing pressure distribution and immersion of the ischial tuberosities. This design delivers a high degree of protection to users who are at High Risk of developing pressure ulcers. This slim cushion format is ideal for both chair and wheelchair use.

Features and Options



Ingress Resistant Cover

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle.

Available in 23 sizes

The *Invacare* **Matrix Flo-tech Image** is available in a range of 23 sizes, from 360 x 330 to 460 x 410 mm. (See size chart for details).

Integral Sag Compensator option

Option of sag compensator to accommodate the sag in wheelchair seats, and ensure the cushion is correctly presented.

Ultra-fresh™ added to all cushions for antimicrobial and odour protection

Antimicrobial treatment that has been added to all **Flo-tech** foam to help make it resistant to the growth of odour causing bacteria. Helps provide longer lasting freshness.

Slimline contoured foam and gel combination cushion






Some users prefer a slimline cushion, finding it more compatible with their daily living activities. The *Invacare* **Matrix Flo-tech Image** combines the benefit of a contoured foam upper surface with a large fluid sac covering the ischial and sacral area, offering increased support and reducing the peak pressures.

Contoured seating surface with gel sac

Offers improved stability and postural support; moulded foam castellations and fluid combine effectively to reduce shear and friction forces on tissue.

Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please see your local *Invacare* website.

Various sizes	Height	Total product weight	Max. user weight ¹	Max. cleaning temperature
				
A range of 50 sizes	78 mm	1.82 kg	140 kg	80°

Manufactured to European Standards: BS EN 1021-1 Cig, BS EN 1021-2 Match and BS 7176 (Crib 5)

1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.

Colours

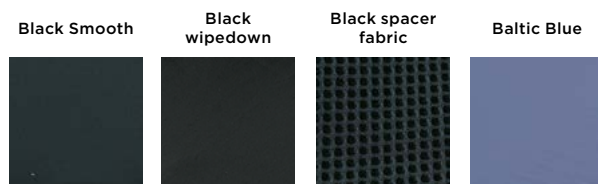

QR Code

Invacare Ltd
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel +44 1 656 776 222
E-mail: uk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS, County Dublin - Ireland
Tel. +353 1 8107084
Email: Ireland@invacare.com



Yes, you can.®

