



Invacare[®] | **Matrix[®] Flo-tech[™] Image**

Slimline moulded foam cushion with fluid sac

matrix[®]
SEATING SERIES

INVACARE
Yes, you can.

The *Invacare Matrix Flo-tech Image* cushion uses a combination of a slim, shaped foam cushion base and a fluid sac covering the sacral area, increasing pressure distribution and immersion of the ischial tuberosities. This design delivers a high degree of protection for people who are at High Risk of developing pressure ulcers. This slim cushion format is ideal for both chair and wheelchair use.

Features and Options



Ingress resistant cover

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle for easy transportation.



Available in 50 sizes

The *Invacare Matrix Flo-tech Image* is available in a range of 50 sizes, from 14-20" width and 13-21" depths. (See prescription chart for details).



Integral sag compensator option

Option to compensate the sag in wheelchair seats, and ensure the cushion is correctly presented.



Ultra-fresh™ added to all cushions for antimicrobial and odour protection

Antimicrobial treatment has been added to all **Flo-tech** foam to help make it resistant to the growth of odour causing bacteria and helps provide longer lasting freshness.

Slimline contoured foam and gel combination cushion

Some individuals prefer a slimline cushion, finding it more convenient with their daily living activities. The *Invacare Matrix Flo-tech Image* combines the benefit of a contoured foam upper surface with a large fluid sac covering the ischial and sacral area, offering increased support and reducing the peak pressures.

Contoured seating surface with gel sac

Offers improved stability and postural support; moulded foam castellations and fluid combine effectively to reduce shear and friction forces on tissue.

Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please visit www.invacare.co.uk

Various sizes	Height	Total product weight	Max. user weight ¹	Max. cleaning temperature	Warranty	Manufactured to European Standards: BS EN 1021-1 Cig, BS EN 1021-2 Match and BS 7176 (Crib 5)
 A range of 50 sizes	 78 mm	 1.82 kg	 140 kg (22 stone)	 80°	 3 years	

1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.

Colours

Black Smooth

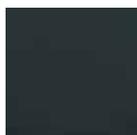

Black wipedown



Black spacer fabric



Baltic Blue

