



Invacare[®] | **Matrix[®] Flo-tech[™] Lite**

Slimline contoured cushion affording exceptional comfort and improved posture

matrix[®]
SEATING SERIES

INVACARE
Yes, you can.

The *Invacare **Matrx Flo-tech Lite*** cushion is a slimline contoured cushion ideally suited for people at low risk of developing pressure ulcers and for comfort purposes only. The contoured seating surface offers improved stability and postural support. With sizes starting from 25 x 25 cm (10 inches), the ***Matrx Flo-tech Lite*** cushion is particularly suited to paediatric clients.

Features and Options



Ingress resistant cover

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle for easy transportation.



Available in 63 sizes

The *Invacare **Matrx Flo-tech Lite*** is available in a range of 63 sizes, from 10-20" width and 10-21" depths. (See prescription chart for details)



Integral sag rigidizer option

Option to compensate the sag in wheelchair seats, and ensure the cushion is correctly presented.



Ultra-fresh™ added to all cushions for antimicrobial and odour protection

Antimicrobial treatment has been added to all ***Flo-tech*** foam to help make it resistant to the growth of odour causing bacteria helps provide longer lasting freshness.

Slimline foam base made from high quality resilient, moulded foam

The high quality resilient foam is manufactured using a unique formula which offers high levels of pressure redistribution and improves product longevity. Moulding allows the ***Matrx Flo-tech Lite*** cushion to provide a 3D form, which cannot be emulated in cut foam. The pommel at the front of the cushion ensures femurs are in a midline position for good posture and allows air to circulate.

Contoured seating surface

Offers improved stability and postural support. Also delivers a greater degree of weight distribution compared to that of a flat surface. The castellations of the cushion allow slight movement which in turn reduces shear and friction forces.

Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please visit www.invacare.co.uk

	Various sizes	Height	Total product weight	Max. user weight ¹	Max. cleaning temperature	Warranty	
							Manufactured to European Standards: BS EN 1021-1 Cig, BS EN 1021-2 Match and BS 7176 (Crib 5)
Matrx Flo-tech Lite	A range of 63 sizes	78 mm	1.5 kg	127 kg (20 stone)	80°	3 years	


1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.

Colours

Black Smooth

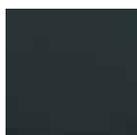

Black wipedown

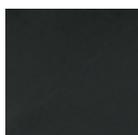

Black spacer fabric



Baltic Blue

