



Invacare®

Matrix® Elan Headrests

Highly-adjustable support and comfort

matrix
SEATING SERIES

INVACARE
Yes, you can.


The elegant *Invacare Matrix Elan* takes our seating range to the next level with its compact, lightweight and adjustable design. The headrest offers superb support and comfort and compliments a full range of paediatric through to adult applications. Key features include:

- Highly adjustable hardware
- Tool-less swingback
- Range of pad options



Highly adjustable hardware

- The innovative hardware on the Elan headrest has six points of adjustment to provide infinite positioning opportunities, including a swing away function, to assist with transfers
- The headrests can be set up quickly and effectively with a single tool while the multi-disk locking feature provides superior holding strength and reliability
- The tool-less adjustment levers make it easy to change the position of the headrest to meet individual's needs
- There are two post options available, the standard adult (380mm) or the mini (280 mm), which can be retro-fitted to Matrix receivers and are compatible with a number of other headrest receivers on the market

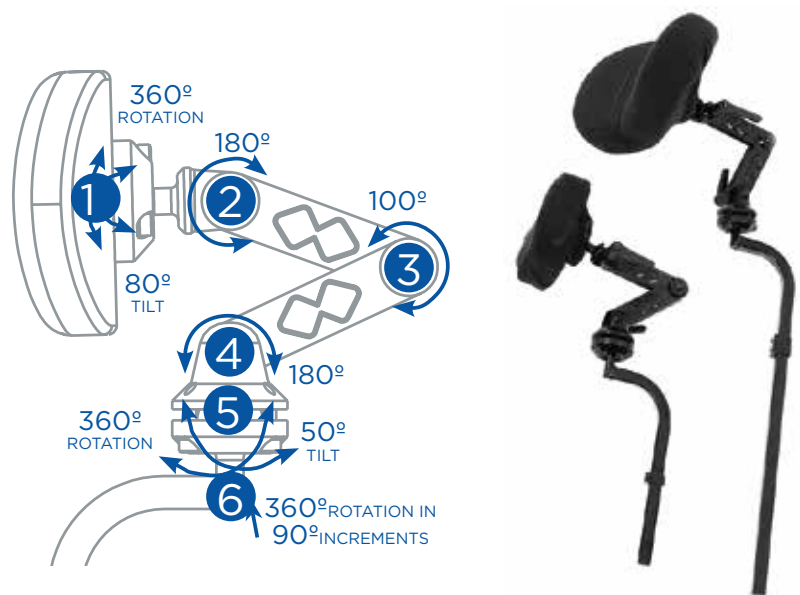

Range of pad options



3 Standard pad options

Ellipse-shaped headrest pads for posterior and lateral support. Can be moulded to fit the shape of the person's head. Available in 3 sizes: 150 x 90 mm, 250 x 130 mm or 360 x 130 mm.



2 Occipital pad options

Provides contoured support to the occipital and sub-occipital areas of the head. Available in two sizes: 230 x 180 mm or 310 x 200 mm.



4-point support

Provides 4 points of support to the occipital and temporal areas of the head for optimum positioning while carefully constructed to fit around the ears. Available in size 280 x 260 mm.

Features

2 tailored cover options

The standard pad cover (Startex fabric) is breathable and sweatresistant. The infection control cover (reverse Startex) is the easy clean version. Both cover options are easy to remove thanks to the elongated dual zip.



Dual density foam

The high density moulded foam provides immersion, comfort, protection and prevents bottoming out. The *Invacare Matrix* pads are 45% lighter than most other pads on the market and surround the edges of the backing plates for extra protection.



Tool-less adjustments

- 1 Adjustment levers allow easy position changes. Fixed screws provided if set position desired.
- 2 Attachment knob has a slim, ergonomic design to allow easy headrest removal or height adjustment.



Interface insert

Two adapters are available which allow the *Matrix Elan* headrest to be fitted to *Invacare Flex 3* solid or tension-adjustable backs.



Multi-disk locking feature

The innovative locking feature provides superior holding strength and reliability.



Adjustable aluminium plates

All pad options can be hand-shaped to fit the needs of the person.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please visit www.invacare.co.uk

Hardware Options

Elan Standard Hardware	380 mm Vertical Post
Elan Mini Hardware	280 mm Vertical Post



matrix[®]
SEATING SERIES

Pad Options

150 mm Ellipse Pad	150 x 90 mm
250 mm Ellipse Pad	250 x 130 mm
360 mm Ellipse Pad	360 x 130 mm
230 mm Occipital Pad	230 x 180 mm
310 mm Occipital Pad	310 x 200 mm
Standard 4 point pad	280 x 260 mm
Small 4 Point pad	230 x 180 mm




Family of seating products

Invacare offer a complete rehabilitation solution: cushion, back support, belt/harness and headrest.



Cover options

Startex



Infection control

