

**matrx**<sup>®</sup>  
SEATING SERIES

# Mini Seating Series



Transforming kids seating



The *Invacare* **Matrx Mini Back** and **Mini PS Cushions** are designed to provide exceptional support with maximum freedom of movement. Every element of the Mini Series has been carefully designed and constructed, keeping paediatric users at the forefront at every stage. Key features include:

- Lightweight, contoured and durable backrest
- Highly resilient contoured foam cushion



#### Lightweight, contoured and durable backrest

The *Invacare* **Matrx Mini Back** provides comfort, stability and improved user positioning. Its Mini Set hardware allows it to be installed and adjusted in minutes while its quick release mechanism makes the lightweight backrest highly portable. The *Invacare* **Matrx Mini Back** offers the possibility to adjust for growth and the adjustable mounting pins allow it to fit to chairs with the same or 25 mm wider back.



#### Highly resilient contoured foam cushion

The *Invacare* **Matrx Mini PS** is an anatomically contoured, high quality foam cushion with a waffled ischial relief section for improved positioning and comfort. The moulded foam provides a three-dimensional form and offers improved stability and support to the user. Moulded foam is also highly durable, increasing the product's longevity.



## Features and Options



### Cover

Back - Water resistant and breathable with attractive Kid\*ab\*ra print.

Cushion - Complete with zippered Startex inner cover and Kid\*ab\*ra print outer cover. It is easy to remove to allow modifications/laundrying.



### Mini Set Back Hardware

Allows 100 mm of depth adjustment to accommodate growth and can recline up to 30° and rotate. The hardware fits 25 mm, 22 mm & 19 mm diameter back canes.



### Mini Swing-away Laterals

Height, width and depth adjustable lateral support for added stability and positioning. Latch release mechanism allows pads to be rotated out of the way for side transfers.

Pad size W110 mm x H82 mm.



### Headrest Clamp

Mini backs are pre-drilled to allow attachment of headrests.



### Mini Elan Headrest

Elan hardware is infinitely adjustable through 6 adjustment points. A full range of Elan positioning pads are available - see Elan headrest brochure for more details.

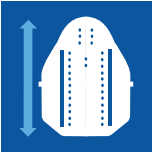
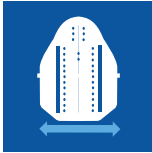
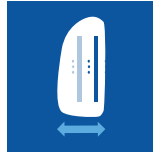











### Chest Harness Interface Kit

Allows attachment of chest support straps to mini backs.

## Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Height	Width	Contour depth	Depth adjustment	Angle adjustment	Total product weight	Max. user weight
Back	 180, 260, 310, 360 mm	 260, 280, 310, 330, 360 mm	 80 mm	 100 mm	 30° range +/- 15°	 400 g and up	 68 kg

	Width	Depth	Height	Total product weight	Max. user weight
Cushion	 260, 280, 310, 330, 360 mm	 260, 310, 360, 410 mm	 50 mm	 150 g and up	 68 kg

## Cushion size availability table

Width (mm) x Depth (mm)	Width (mm) x Depth (mm)
260 x 260 mm	330 x 310 mm
260 x 310 mm	330 x 360 mm
260 x 360 mm	330 x 410 mm
280 x 310 mm	360 x 310 mm
280 x 360 mm	360 x 360 mm
280 x 410 mm	360 x 410 mm
310 x 310 mm	
310 x 360 mm	
310 x 410 mm	



## Cover options

Kid\*ab\*ra print



QR Code

Invacare Ltd  
Pencoed Technology Park - Pencoed  
Bridgend CF35 5AQ - United Kingdom  
Tel +44 1 656 776 222  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS, County Dublin - Ireland  
Tel. +353 1 8107084  
Email: [Ireland@invacare.com](mailto:Ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

