

matrix[®]
SEATING SERIES

Mini Seating Series



Transforming kids seating



The *Invacare* **Matrx Mini Back** and **Mini PS Cushions** are designed to provide exceptional support with maximum freedom of movement. Every element of the Mini Series has been carefully designed and constructed, keeping paediatrics at the forefront at every stage. Key features include:

- Lightweight, contoured and durable backrest
- Highly resilient contoured foam cushion
- Highly-adjustable **Matrx Mini Elan** headrest



Lightweight, contoured and durable backrest

The *Invacare* **Matrx Mini Back** provides comfort, stability and improved positioning. Its Mini Set hardware can be installed and adjusted in minutes while its quick release mechanism makes the lightweight backrest highly portable.

The *Invacare* **Matrx Mini Back** offers the possibility to adjust for growth and the adjustable mounting pins allow it to fit to chairs with the same or 25 mm wider back.



Highly resilient contoured foam cushion

The *Invacare* **Matrx Mini PS** is an anatomically contoured, high quality foam cushion with a waffled ischial relief section for improved positioning and comfort. The moulded foam provides a three-dimensional form and offers improved stability and support to the individual. Moulded foam is also highly durable, increasing the product's longevity.



Features and Options



Cover

Back - Water resistant and breathable with attractive Kid*ab*ra print.

Cushion - Complete with zippered Startex inner cover and Kid*ab*ra print outer cover. It is easy to remove to allow modifications/laundrying.



Mini Set Back Hardware

Allows 100 mm of depth adjustment to accommodate growth and can recline up to 30° and rotate. The hardware fits 25 mm, 22 mm & 19 mm diameter back canes.



Mini Swing-away Laterals

Height, width and depth adjustable lateral support for added stability and positioning. Latch release mechanism allows pads to be rotated out of the way for side transfers. Pad size W110 mm x H82 mm.



Headrest Clamp

Mini backs are pre-drilled to allow easy attachment of headrests.



Mini Elan Headrest

Elan hardware is infinitely adjustable through 6 adjustment points with a full range of Elan positioning pads - see the Elan headrest brochure for further details.




Chest Harness Interface Kit

Allows attachment of chest support straps to mini backs.

Cover options

Kid*ab*ra print

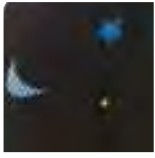

Startex inner cover with funky print on outer cover

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please visit www.invacare.co.uk

	Height	Width	Contour depth	Depth adjustment	Angle adjustment	Total product weight	Max. user weight
Back	180, 260, 310, 360 mm	260, 280, 310, 330, 360 mm	80 mm	100 mm	30° range +/- 15°	400 g and up	68 kg

	Width	Depth	Height	Total product weight	Max. user weight
Cushion	260, 280, 310, 330, 360 mm	260, 310, 360, 410 mm	50 mm	150 g and up	68 kg (10.7 stone)

Back

		WIDTH				
		26	28	31	33	36
DEPTH	cm					
	inch	10	11	12	13	14
	18	07				
	26	10				
	31	12				
36	14					

Cushion

		WIDTH				
		26	28	31	33	36
DEPTH	cm					
	inch	10	11	12	13	14
	26	10				
	31	12				
	36	14				
41	16					

