


Essential Touch

Double-sided castellated foam mattress

A practical support surface offering
comfort, and pressure redistribution.

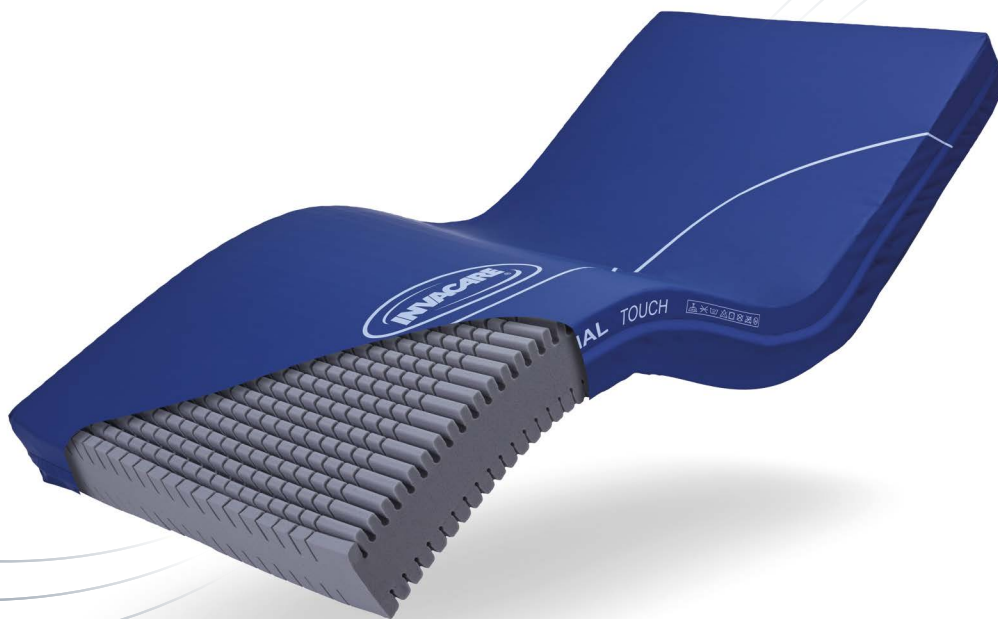

A cost-effective solution

For pressure injury prevention



Yes, you can.®

Affordable quality mattress



Pressure Redistribution

- › Keyhole cut cells move independently to distribute weight evenly.
- › Castellated foam with independent cells offer greater immersion and envelopment.
- › Larger outer castellations and smaller inner castellations where pressure is more prominent.



Shear & Friction Reduction

- › Castellations help reduce shear when repositioning on the bed.
- › Multi-stretch, water-resistant, polyurethane outer cover further helps to reduce shear and friction.
- › Suitable for multiple care environments.



Comfort

- › Castellations provide greater comfort and immersion.
- › Double-sided castellated foam allows the mattress to be turned or rotated inside the cover for increased product longevity.
- › Cover with two reinforced handles and zipped around three sides for easy handling and removal for laundering.

Technical data ▼

For more information about this product, including the user manual, please visit our website: www.invacare.co.uk

Sizes available ▼

		830 mm	850 mm	880 mm	900 mm	1030 mm	1180 mm
150 mm height	2000 mm	●	●	●	–	●	●
	2050 mm	–	–	●	–	–	–
	2200 mm	–	–	–	●	–	–



Max. User Weight ▼

150 kg



Max. Cleaning Temperature ▼

95 °C

Manufactured to: BS EN 597-1 Cig, BS EN 597-2 Match and BS 7177 (ignition Source 5)

Essential nursing care is pivotal in pressure injury management. This mattress may positively contribute to the outcome of a pressure injury prevention plan that includes education, clinical judgment and action planning based on vulnerability, which are fundamental factors in the prevention of pressure injuries. A range of assessment scales can be used as a formal method of assessing risk from pressure injury development, and should be used in conjunction with an informal assessment (informed nursing judgement) as part of a comprehensive and individualised plan.

T: +44 1656 776222
uk@invacare.com
www.invacare.co.uk
www.invacare.ie

Invacare UK & Ireland
Unit 4 Pencoed Technology
Park
Bridgend CF35 5AQ



Yes, you can.®



© 2024 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. Essential Touch - UK - 11/2024