



Invacare®

Action[®] 3_{NG} Light

Outstanding performance



Introducing the all new **Action³_{NG} Light**. Elevating the **Action** range even further with new cutting edge designs that enable the chair to smoothly glide - delivering effortless propulsion. **Action³_{NG} Light** is built upon a lightweight aluminium frame that is designed to compliment daily activities at home and out and about - mobility made easy.



- ✓ From as light as 13.0kg
- ✓ Quality performance
- ✓ Sleek black components
- ✓ Lightweight rear wheels
- ✓ Improved positioning and comfort

Action³_{NG} Light



New innovative lightweight castor housing greatly reduces the overall weight of the chair with its sleek aluminium design, angle adjustable from 0°, 3° to 6°.

Setting the adjustment is easier than ever before with a clear two part plate process enabling seat angles to be adjusted for individual positioning.

Transport weight from - 7.1kg



The **Action³_{NG} Light** is designed with both the individual and carer/attendant in mind - effortless to lift and easy to transport. The lightweight aluminium frame allows for easier propulsion and energy conservation. Light to lift making transportation easy with removable features and a compact frame.

Features and options



A new angle adjustable backrest provides precise positioning from -15° to $+15^{\circ}$, which can be carried out in one smooth process whilst the individual remains in the chair.



Folding is quick and simple with an innovative pull cord allowing the backrest to fold directly on top of the seat.



Newly designed black mesh bimaterial back upholstery delivers increased comfort and postural support. Shaped to overlap the back canes and features a new removable lumbar strap for optimal positioning.



20 degree centre of gravity positions



Black front and rear wheels



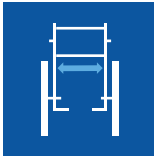






Bodypoint 2 point positioning belt


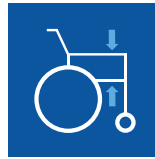
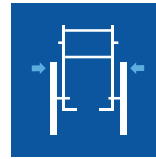











A choice of cushions individually designed to relieve pressure and provide comfort and support.

Technical data

For more comprehensive information about this product, including the product's user manual, please visit www.invacare.co.uk

	Seat width	Seat depth	Seat angle	Seat height	Backrest angle	Backrest height	Front frame angle
<i>Invacare</i> Action³_{NG} Light	 380 / 405 / 430 / 455 / 480 / 505 mm	 400 / 450 / 500 mm	 0° / 3° / 6°	 385 / 410 / 435 / 460 / 485 / 510 mm	 0° / 7° / -15° - +15°	 325 - 510 mm	 80° / 0 - 80° / 90°

	Length of footrests	Armrest height	Driving unit width	Total width, folded	Total length inc. footrests	Total length without footrests	Total height
<i>Invacare</i> Action³_{NG} Light	 200 - 400 / 330 - 460 / 440 - 510 mm	 200 - 280 / 240 / 170 - 300 / 210 - 350 mm	 Seat width + 190 mm	 295 mm	 980 - 1105 mm	 835 - 960 mm	 820 - 1050 mm

	Max. safe slope	Radius	Max. user weight	Total weight from	Mass of the heaviest part, from
<i>Invacare</i> Action³_{NG} Light	 9°	 870 mm	 125 kg	 13.0 kg	 7.1 kg



Frame colours

Hypnotic Purple



Matt Black



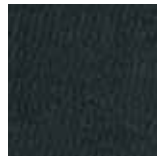
Blue Pop




Matt Red




Black Nylon



Pearl White



Slate Grey



Electric Green



Anthracite Black



Upholstery



Please notice colours may vary slightly from those displayed above.

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

Trusted Solutions,
Passionate People

