



rea

Confidence in posture solutions

rea | **Azalea[®] Max**

A tailor-made solution for heavier individuals

INVACARE
Yes, you can.

The **Rea Azalea Max** is designed to meet the specific needs of people with extra weight. For effortless positioning we recommend power tilt and backrest recline. The **Rea Azalea Max** is a niche product which complements the existing **Azalea** family from Invacare.

Continued Compatibility

The **Rea Azalea Max** includes features such as tilt-in-space with weight-shifting and DSS (Dual Stability System) mechanism from the **Azalea** range. Stability in this wheelchair originates from body length compensation in the back and legrests, beneficial for those with limited movement who have difficulty changing position.

Features and Options



Two-way backrest

The lower and upper parts of the shaped backrest are reversible, giving the option of more space in the seat or chest area.



Footplate converter

Divided into three parts, the footplate comes with a removable mid section, to facilitate transfers.



Power tilt and recline

With the hand control the person or the attendant can easily change seating position.



High ergonomic push handle

The ergonomically shaped push handles are positioned to optimise the pushing for the attendant. An alternate version with higher handles offers a different position for the hands and is available as an option.

Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please visit www.invacare.co.uk

	Seat width	Seat depth	Seat height	Backrest height	Armrest height	Length of legrests	Seat angle	Backrest angle
Rea Azalea Max	550; 610 + 50 mm	530 - 600 mm	400 / 450 mm	500 + 80 mm	320 - 420 mm	330 - 500 mm	1° - 20°	0°-30°
Rea Azalea Max Assist	550; 610 + 100 mm	530 - 600 mm	400 / 450 mm	500 + 80 mm	320 - 420 mm	330 - 500 mm	1° - 20°	0°-30°


	Driving unit width	Total height	Total length inc. Footrests	Total weight	Max. User weight	Transport weight	Maximum slope for use of parking brake
Rea Azalea Max	SW + 250 mm	1020 - 1240 mm	1050 - 1120 mm	54 kg	max. 180 kg	33 kg	7°
Rea Azalea Max Assist	SW + 220 mm	1020 - 1240 mm	1050 - 1120 mm	54 kg 1	max. 180 kg	37 kg	7°

1. Transport weight increases by 4,2 kg if Power tilt & recline is selected.

Frame colour

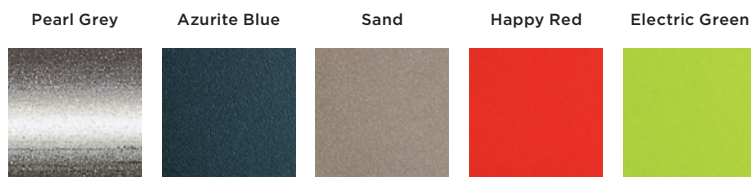

Upholstery

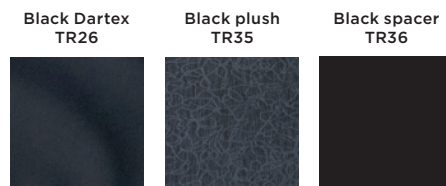

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

